Exporting QRDA I Files from DrChrono

07/08/2024 7:44 pm EDT


You may need to export a QRDA I file to provide specific evidence for a particular patient. Or you may need to export your Clinical Quality Measure data for all the measures available in DrChrono or specific ones.

This article will show you how to export QRDA I files by measure(s) and by patient.


Exporting a QRDA by Measure(s)

Go to Clinical > Clinical Quality Measures.


On the main page, select the reporting year or date range for your export.

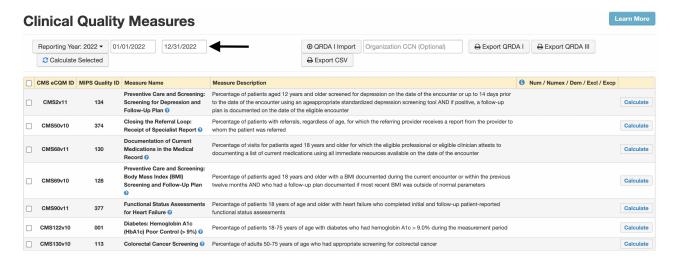

Next, select the measure or measures you would like to include in your export. Click **Calculate Selected** or **Calculate** next to the individual measure to view the data associated with the measure.

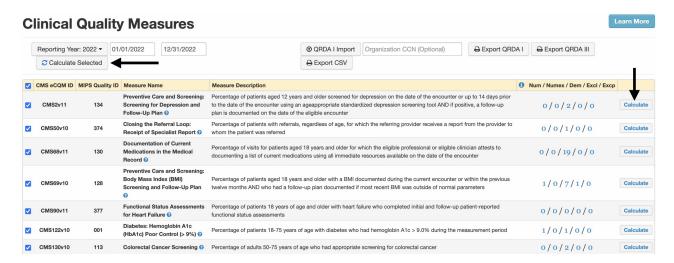


Click Export QRDA I to export the QRDA file.

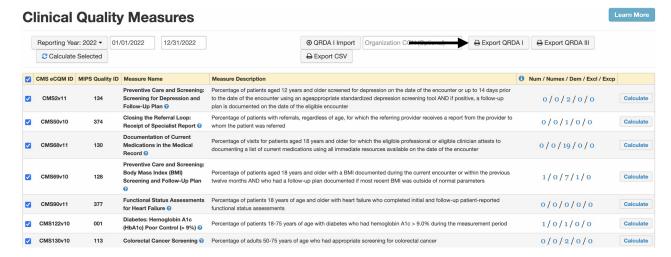

The file will export to the Message Center (




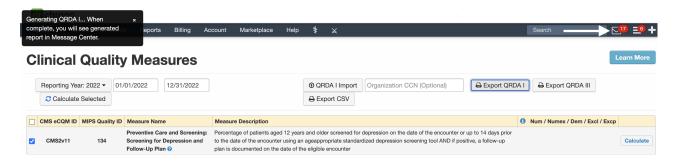

In the Message Center, you can select download next to the file.




You can also open an individual message and download it from there.

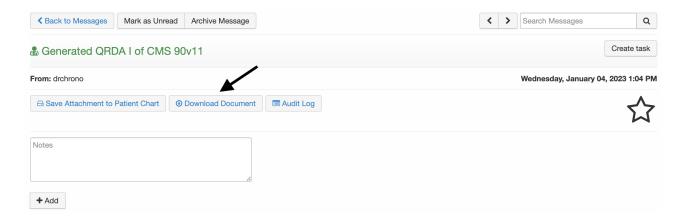

Exporting a QRDA I for an Individual Patient

The QRDA I export process for an individual patient is straightforward.

1. Go to Patients > Patient List and select the patient or search for the patient you need.

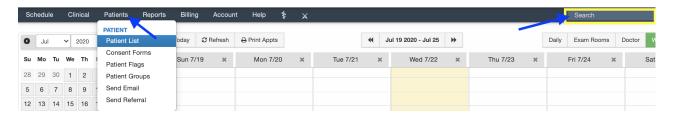

- 2. Select CQMs from the left navigation panel, Next, click Download QRDA I.
- 3. Save the QRDA I XML file based on your computer's settings.

